


Askina® – Quality has a name

Dressings for veterinary medicine

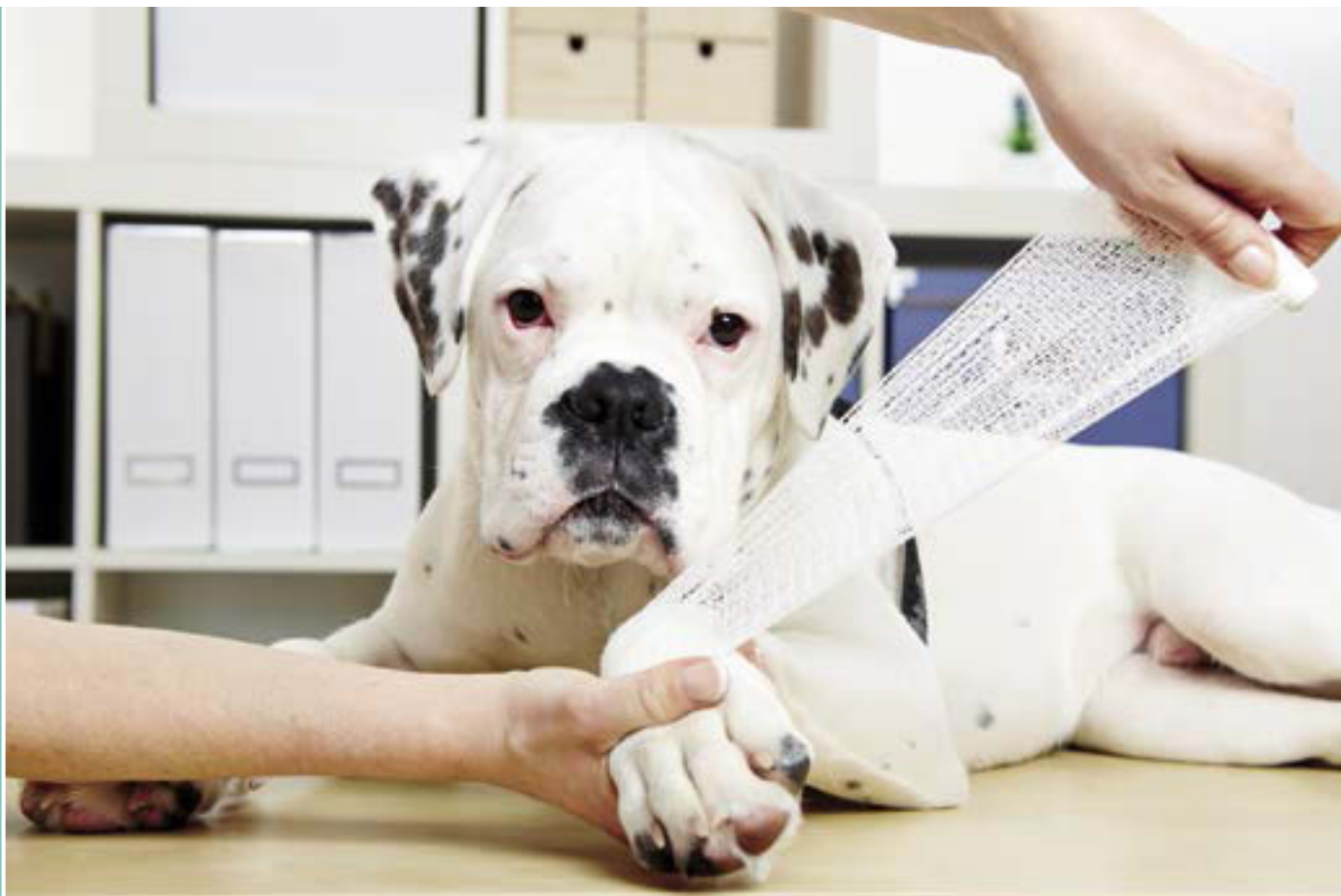

B. Braun Vet Care

Only a clean wound can heal

B. Braun offers an integrated concept for modern wound care, as well as a comprehensive range for traditional wound care.



Traditional wound care

Adhesive dressing material

Askina® Plast E	5
Askina® Vet	6

Cohesive dressing material

Askina® Protect	7
Askina® Fix E	8

Cotton wool padding

Askina® Soft Pad	9
------------------	---

Tricot tube

Askina® TS-Bandage	10
--------------------	----

Cast

Askina® Cast	11
--------------	----

Compression bandages

Askina® Gaze compresses	12
-------------------------	----

Fixation bandages

Askina® Silk	13
Askina® Film	13

Modern wound care

Wound irrigation solution and antiseptic

ProntoVet® Solution / Gel	14
---------------------------	----

Wound paste

Askina® Calgitrol®	17
--------------------	----

Wound dressings

Askina® Sorb	18
Askina® SilNet	19
Askina® DresSil / DresSil Border	20
Askina® Foam	25

Dependable, strong and functional.

Askina® products for traditional wound care



B | BRAUN
SHARING EXPERTISE

Wound care and wound preparation

B. Braun understands the need for proper wound care. We offer a product range that helps to achieve these goals – specifically for veterinary medicine.

	Adhesive		Cohesive		Padding	Underlay	Fixation
	Askina® Plast E	Askina® Vet	Askina® Protect	Askina® Fix E	Askina® Soft Pad	Askina® TS-Bandage	Askina® Cast
Fixating			✓	✓		✓	
Light cushioning	✓	✓	✓				
Firm cushioning	✓	✓	✓			✓	✓
Immobilizing	✓	✓					✓
Cotton wool padding					✓		
Self adherent	✓	✓	✓	✓			✓
Elastic	✓	✓	✓	✓		✓	
Soft				✓	✓	✓	
Breathable				✓	✓		
Easily moldable to small body parts				✓	✓		
Water-repellant	✓	✓	✓				✓
Easy to tear			✓		✓		
High cotton content	✓			✓			

Askina® Plast E

Premium adhesive bandage



Adhesive dressing material

Askina® Plast E is a premium quality product in the range of adhesive bandages, made from 100% cotton and air- and vapor permeable due to its zinc oxide coating.

- High wear comfort due to natural material properties
- Strong, green plastic core
- Individually packaged
- Pure white – for clean, functional appearance

Practical tips on application

- Ideal for bandages that require a tight hold

Indications

- Specially designed for equestrian sports
- Sling bandages
- Firm compression bandages

Art. No.	Description	Width	Length (stretched)	Pcs/Pack
9058005DE	Askina® Plast E	5 cm	4.5 m	12 pcs
9058006DE	Askina® Plast E	7.5 cm	4.5 m	12 pcs
9058007DE	Askina® Plast E	10 cm	4.5 m	12 pcs

Askina® Vet

Standard adhesive bandage



Askina® Vet is a standard adhesive bandage made from a polyurethane-cotton fabric that allows for high elasticity.

- Strong, green plastic core
- Individually packaged

Practical tips on application

- Economical and strong


Indications

- Sling bandages
- Firm compression bandages

Article number	Description	Width	Length (stretched)	Pcs/Pack
9051075DE	Askina® Vet	7.5 cm	4.5 m	12 pcs

Askina® Protect

Universal cohesive bandage

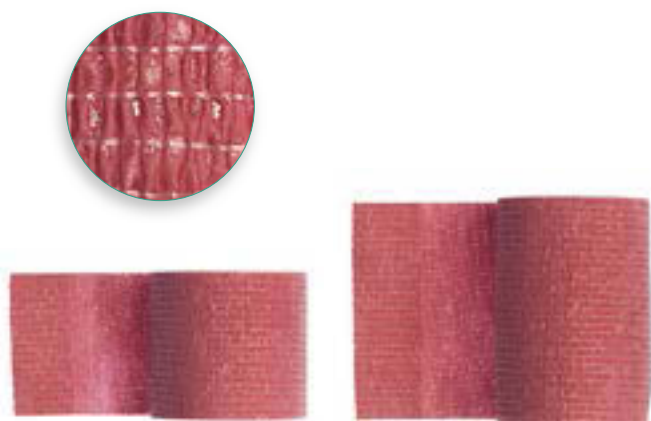

Cohesive dressing material

Askina® Protect is a multifunctional synthetic bandage.

- Very good cohesive properties
- Light and very comfortable for patients
- Individually packaged
- Tightly woven fabric structure for an optimal hold
- Water-repellant

Practical tips on application

- Tightly woven fabric structure – less need of material

Indications

- Fixation for dressings and bandages
- Fixation for padding
- Light compression bandages

Article number	Description	Width	Length (stretched)	Pcs/Pack
9057050DE	Askina® Protect	5 cm	4.5 m	12 pcs
9075075DE	Askina® Protect	7.6 cm	4.5 m	12 pcs

Askina® Fix E

Premium cohesive fixation fleece



Askina® Fix E is a premium-quality cohesive fixation fleece with high cotton content. Askina® Fix E allows for an outstanding heat- and fluid exchange which facilitates the healing process.

- Excellent cohesive properties
- Sterilizable
- Individually packaged

Practical tips on application

- The time-saving alternative to gauze bandages


Indications

- Fixation of wound bandages, IV catheters and wound dressings
- For use with Robert Jones bandages

Article number	Description	Width	Length (stretched)	Pcs/Pack
9059145DE	Askina® Fix E	4 cm	4.0 m	12 pcs
9059148DE	Askina® Fix E	10 cm	4.0 m	12 pcs

Askina® Soft Pad

Cotton wool padding

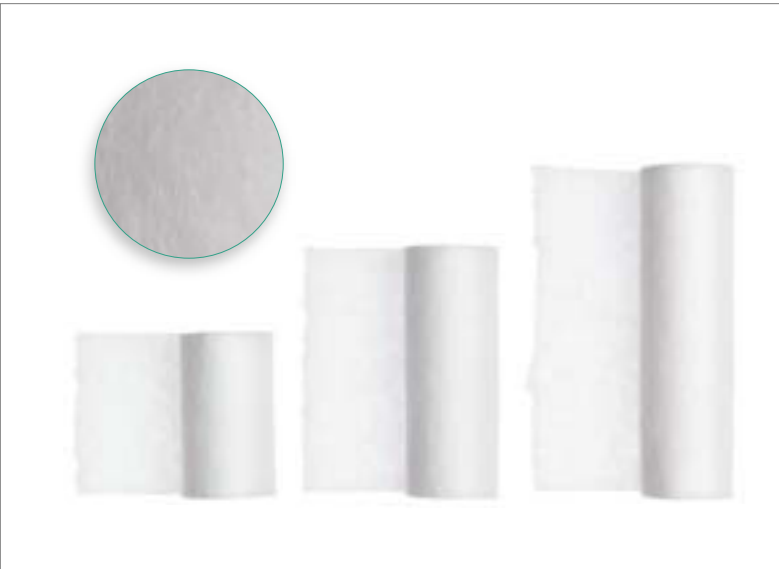

Cotton wool padding

Askina® Soft Pad is a light, breathable cotton wool padding ideal for the padding of dressings.

- Breathable and permeable to exudate
- Individually packed in paper banderole

Practical tips on application

- Exceptional stability and easy to tear

Indications

- For primary padding
- Underlay for plaster-, compression- and rigid cast bandages as well as splints

Article number	Description	Width	Length (stretched)	Pcs/Pack
9058118DE	Askina® Soft Pad	10 cm	3 m	12 pcs
9058126DE	Askina® Soft Pad	15 cm	3 m	12 pcs
9058134DE	Askina® Soft Pad	20 cm	3 m	6 pcs

Askina® TS-Bandage

Tricot tube



Askina® TS-Bandage is an easy to apply tricot tube that saves time.

- Skin-friendly underlay bandage made of 100 % cotton
- Easy to apply
- Can be sterilized if necessary
- Elastic
- Run-resistant

Practical tips on application

- Fast and cost-saving

Indications

- As a crease-free underlay for plaster dressings / synthetic protective bandages
- Fixation of wound dressings
- Trunk- and head dressings

Article number	Description	Size	Width	Length (stretched)	Pcs/Pack
9043314	Askina® TS-Bandage	1	2.5 cm	20 m	20 pcs
9043322	Askina® TS-Bandage	2	4.0 cm	20 m	20 pcs
9043330	Askina® TS-Bandage	3	6.0 cm	20 m	10 pcs
9043349	Askina® TS-Bandage	4	8.0 cm	20 m	10 pcs
9043357	Askina® TS-Bandage	5	10.0 cm	20 m	10 pcs
9043365	Askina® TS-Bandage	6	12.0 cm	20 m	10 pcs
9043373	Askina® TS-Bandage	7	16.0 cm	10 m	10 pcs

Askina® Cast

The synthetic support bandage



Cast

Askina® Cast is a synthetic support bandage for stable support and immobilization of limbs.

- Easy and hassle-free to unroll
- Easy modelling
- Radiolucent
- Water-resistant (can be washed off)
- Fully hardened after approximately 15 minutes
- Fully loadable after approximately 30 minutes
- Individually packaged

Practical tips on application

- No production of heat

Indications

- Treatment of fractures
- Post-operative immobilization

Article number	Description	Width	Length	Pcs/Pack
9073108	Askina® Cast	5 cm	3.6 m	10 pcs
9073116	Askina® Cast	7.6 cm	3.6 m	10 pcs
9073124	Askina® Cast	10.1 cm	3.6 m	10 pcs

Askina® Gauze Compresses

Compresses for all indications



- Gauze compresses made of absorbent gauze with folded cut edges
- Bleached 100 % without chlorine
- 8-layered and 12-layered
- Available sterile and non-sterile packed
- 100 % cotton

Article number	Description	Size	Pcs/Pack
9031103N	Askina® Gauze Compresses 8-layered, sterile	5 x 5 cm	5 x 2 pcs
9031111N	Askina® Gauze Compresses 8-layered, sterile	7.5 x 7.5 cm	5 x 2 pcs
9031120N	Askina® Gauze Compresses 8-layered, sterile	10 x 10 cm	5 x 2 pcs
9031138N	Askina® Gauze Compresses 8-layered, sterile	10 x 20 cm	5 x 2 pcs
9031200N	Askina® Gauze Compresses 8-layered, sterile	5 x 5 cm	25 x 2 pcs
9031219N	Askina® Gauze Compresses 8-layered, sterile	7.5 x 7.5 cm	25 x 2 pcs
9031227N	Askina® Gauze Compresses 8-layered, sterile	10 x 10 cm	25 x 2 pcs
9031235N	Askina® Gauze Compresses 8-layered, sterile	10 x 20 cm	25 x 2 pcs
9031308	Askina® Gauze Compresses 8-layered, nonsterile	5 x 5 cm	100 pcs
9031316	Askina® Gauze Compresses 8-layered, nonsterile	7.5 x 7.5 cm	100 pcs
9031324	Askina® Gauze Compresses 8-layered, nonsterile	10 x 10 cm	100 pcs
9033017	Askina® Gauze Compresses 12-layered, nonsterile	7.5 x 7.5 cm	100 pcs
9033025	Askina® Gauze Compresses 12-layered, nonsterile	10 x 10 cm	100 pcs
9033041	Askina® Gauze Compresses 12-layered, nonsterile	10 x 20 cm	100 pcs

Askina® Silk /Film

Bandages for the perfect fixation of dressings



Fixation bandages

Askina® Silk

- Reliable adhesion on different surfaces and optimal skin friendliness
- Modern acrylate adhesive
- Easy to tear
- Removal with little pain

Askina® Film

- Self-adhesive, breathable fleece
- Very moldable and elastic
- Joints can be bandaged without restriction of movement
- Easy fixation due to two-part protective film

Indications

- Fixation of wound dressings of all types
- Fixation of catheters, tubes and puncture sets

Article number	Description	Width	Length	Pcs/Pack
9080112	Askina® Silk	1.25 cm	9.1 m	24 pcs
9080120	Askina® Silk	2.5 cm	9.1 m	12 pcs
9080015	Askina® Silk	1.25 cm	5 m	24 pcs
9080023	Askina® Silk	2.5 cm	5 m	12 pcs
9080031	Askina® Silk	5 cm	5 m	6 pcs
9080228	Askina® Silk	2.5 cm	5 m	1 pcs
9081500	Askina® Film	1.25 cm	9.1 m	24 pcs
9081518	Askina® Film	2.5 cm	9.1 m	12 pcs

Individual, phase-specific and simple.

Askina® products for modern wound care



B | BRAUN
SHARING EXPERTISE

ProntoVet®

The premium quality irrigation solution specifically for veterinary medicine



Wound Irrigation Solution and Antiseptic

Only a clean wound can heal.

The patented ProntoVet® formulation of polihexanide and betaine cleans the wound and gently removes biofilm.

ProntoVet® is non-irritating and pain-free, very well tolerated by skin and mucous membranes, and acts without restriction in the presence of blood and protein.

- Sterile
- Colorless
- Ready for use
- Can be used for 8 weeks after opening*
- Practical cap with sealing pin
- Gel consistency prevents flow off from the target location (ProntoVet® Gel)










* at room temperature

Indications

- Irrigation of wounds and wound cavities of all types
- Gentle to tissue and pain-free removal of crusted dressings
- Moistening of dressings
- Long lasting, moist wound covering (ProntoVet® Gel)

Article number	Description	Size	Pcs/Pack
400600	ProntoVet® woundgel	30 ml	20 pcs
400601	ProntoVet® wound irrigation solution	350 ml	10 pcs

Askina® – Focus products at one glance

	Abrasions	Traumatic wounds / Cut and stab wounds	Surgical wounds	Burns	Thermal and chemical wounds	Infected wounds	
Askina® SilNet Soft silicone backing with safe adhesion	✓	✓	✓	✓	✓		  <ul style="list-style-type: none"> ■ Atraumatic dressing changes ■ Easy to apply ■ Can stay on the wound for several days ■ Hypoallergenic
Askina® Calgitrol® Silver alginate in a tube	✓	✓	✓	✓	✓	✓	 <ul style="list-style-type: none"> ■ Simple and effective application ■ Very good adaptability to wound contours ■ Controlled, steady liberation of silver ions ■ Active against a broad spectrum of microorganisms, including MRSA ■ Easy removal by rinsing
Askina® DresSil Border Foam dressing with soft silicon adhesive and adhesive border	✓	✓	✓	✓	✓		 <ul style="list-style-type: none"> ■ Atraumatic dressing changes ■ High level of comfort for the patient ■ Can stay on the wound for several days ■ Usable under compression bandages ■ Hypoallergenic
Askina® Sorb Conformable alginate dressing / -tamponade	✓	✓	✓	✓	✓		 <ul style="list-style-type: none"> ■ Cleans and protects the wound edges ■ Conformable and removable in one piece ■ Does not leak and maintains a moist wound environment ■ Absorbs excess exudate ■ Suitable for prolonged use
Askina® Foam Highly absorbent foam dressing for heavy exudation	✓	✓				✓	 <ul style="list-style-type: none"> ■ Excellent fluid handling characteristics ■ Good retention capacity ■ Usable under compression bandages
<div>  Low weeping wound – Little exudate </div> <div>  Moderately weeping wound – Moderate exudate </div> <div>  Heavily weeping wound – Heavy exudate </div>							

✓ especially suitable for ...

Askina® SilNet

Soft silicone wound contact layer



- Atraumatic dressing changes
- Easy to use
- Provides very good management of exudate
- Can stay on the wound for several days
- Does not stick to wound

Practical tips on application

Askina® SilNet can be combined with a variety of secondary bandages and dressings.

Indications

- Traumatic wounds, e.g. cuts or abrasion wounds
- Thermal wounds, e.g. burns
- Post-operative wounds
- Pre-damaged skin
- Care of skin grafts and donor sites

Article number	Description	Size	Pcs/Pack
5195710	Askina® SilNet	5 x 7.5 cm	10 pcs
5195703	Askina® SilNet	5 x 7.5 cm	3 pcs
5197510	Askina® SilNet	10 x 7.5 cm	10 pcs
5197503	Askina® SilNet	10 x 7.5 cm	3 pcs
5191810	Askina® SilNet	10 x 18 cm	10 pcs
5191803	Askina® SilNet	10 x 18 cm	3 pcs
5192305	Askina® SilNet	20 x 30 cm	5 pcs

Askina® Calgitrol®

Silver alginate in a tube



- Simple, targeted application
- Controlled silver ion delivery for long lasting antibacterial protection
- Active against a broad spectrum of microorganisms, including MRSA
- No activation needed
- Highly conformable paste
- Maintains a moist wound environment
- Easy to remove (ProntoVet®) from the wound by simple rinsing

Practical tips on application

- Knead or shake tube prior to use
- Removal of Askina® Calgitrol® by thorough irrigation with ProntoVet®

Indications

- Burns
- Infected wounds
- Traumatic wounds

Article number	Description	Size	Pcs/Pack
6241505	Askina® Calgitrol®	15 g	5 pcs
6241510	Askina® Calgitrol®	15 g	10 pcs

Askina® DresSil / DresSil Border

Foam dressing with silicone adhesive



- Optimal exudate management
- Easy to use
- High level of comfort for the patient
- Can stay on the wound for several days
- Usable under compression bandages

Indications

- Surgical wounds
- Traumatic wounds, e.g. cuts and abrasions
- Heavily exuding wounds


Practical tips on application

The wound dressings are coated with a soft silicone adhesive, providing gentle and secure adherence. The silicone adapts to the uneven contours of the skin and forms a larger contact area compared to traditional adhesion bases.

Article number	Description	Size with adhesive border	Size without adhesive border	Pcs/Pack
5295703	Askina® DresSil		5 x 7 cm	3 pcs
5295710	Askina® DresSil		5 x 7 cm	10 pcs
5291003	Askina® DresSil		10 x 10 cm	3 pcs
5291010	Askina® DresSil		10 x 10 cm	10 pcs
5291503	Askina® DresSil		15 x 15 cm	3 pcs
5291510	Askina® DresSil		15 x 15 cm	10 pcs
5291203	Askina® DresSil		10 x 20 cm	3 pcs
5291210	Askina® DresSil		10 x 20 cm	10 pcs
5292005	Askina® DresSil		20 x 20 cm	5 pcs
5396610	Askina® DresSil Border	6 x 6 cm	3 x 3 cm	10 pcs
5397503	Askina® DresSil Border	7.5 x 7.5 cm	4.5 x 4.5 cm	3 pcs
5397510	Askina® DresSil Border	7.5 x 7.5 cm	4.5 x 4.5 cm	10 pcs
5391003	Askina® DresSil Border	10 x 10 cm	7 x 7 cm	3 pcs
5391010	Askina® DresSil Border	10 x 10 cm	7 x 7 cm	10 pcs
5391503	Askina® DresSil Border	15 x 15 cm	12 x 12 cm	3 pcs
5391510	Askina® DresSil Border	15 x 15 cm	12 x 12 cm	10 pcs
5395203	Askina® DresSil Border	15 x 20 cm	12 x 17 cm	3 pcs
5395210	Askina® DresSil Border	15 x 20 cm	12 x 17 cm	10 pcs
5392005	Askina® DresSil Border	20 x 20 cm	17 x 17 cm	5 pcs

Askina® Sorb

Alginate dressing with excellent stability

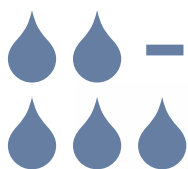

- Flexible and easy to apply
- Conforms to the contours of the wound
- High absorption capacity
- Can be removed in one piece
- Exceptional stability due to linked fiber structure
- Atraumatic dressing changes with little pain
- Woven alginate-carboxymethyl cellulose (CMC) fiber structure

Practical tips on application

In case of dry wound conditions, use in combination with ProntoVet® wound irrigation solution and ProntoVet® wound gel.

Tips on choosing the proper dressing:

- Flat fixation with non-woven sheet dressings and marginal fixation with adhesive tape
- Circular fixation with fixation bandages from the traditional product range

Indications

- Thermal / chemical wounds, burns
- Abrasions
- Traumatic wounds, surgical wounds

Article number	Description	Size	Pcs/Pack
2109S	Askina® Sorb	6 x 6 cm	3 pcs
2115S	Askina® Sorb	6 x 6 cm	10 pcs
2107S	Askina® Sorb	10 x 10 cm	3 pcs
2116S	Askina® Sorb	10 x 10 cm	10 pcs
2102S	Askina® Sorb	15 x 15 cm	10 pcs
2106S	Askina® Sorb Tamponade	2.7 x 34 cm	3 pcs

Standard dressings in veterinary medicine

Attaching an IV catheter

Askina® Silk | Askina® Film

Askina® Fix E



Fixating a catheter to skin with adhesive tape



Fixating a catheter to a leg with fixation bandage

Padding

Askina® Soft Pad

Askina® Soft Pad

Askina® Fix E



Applying padding between toes and joints



Fully wrapping the wound site



Cotton wool with fixation bandage

Askina® Fix E


Askina® Protect



Tightening a fixation bandage



Fixating the dressing with a cohesive bandage

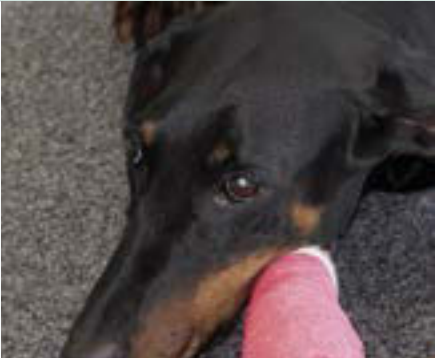

Robert Jones bandage

Askina® Silk | Askina® Film | Askina® Soft Pad

Askina® Soft Pad



Fixating cotton wool on a paw with adhesive tape



Wrapping the entire leg with padding

Askina® Fix E

Askina® Protect



Applying cotton wool padding with fixation bandage



Fixating the dressing with a cohesive bandage

Guide to chronic wound care

Wound description / stage of wound healing					Therapeutic goal							
		Dry	Moist	Wet	Improvement of blood circulation	Removal of necrotic tissue	Removal of necrotic and fibrinous tissue	Wound base free of fibrinous tissue, promotion of granulation	Promotion of granulation	Promotion of granulation and epithelialization	Stabilization of newly formed tissue	Elimination of bacterial load
Black necrotic wound					✓	✓						
Black necrotic and fibrinous wound		✓				✓						
Black necrotic, fibrinous and granulating wound			✓	✓			✓					
Fibrin coatings		✓						✓				
Fibrin coatings and granulation			✓	✓				✓				
Granulation		✓							✓			
			✓	✓					✓			
Granulation and epithelialization		✓								✓		
			✓	✓						✓		
Epithelialization		✓	✓								✓	
Infected wounds		✓										✓
			✓	✓								✓

Procedure / actions	Wound cleansing /- irrigation				Wound filler			Wound dressing			
	Cleansing with wound irrigation solution	Soaking with wound irrigation solution, and mechanical cleaning	Wound cleansing without aggressive mechanical cleaning	Wound irrigation	Wound Gel	Askina® Sorb	Askina® Calgitrol®	Askina® SiINet	Askina® DresSil	Askina® DresSil Border	Askina® Foam
Surgical debridement (then new wound classification), rehydration, promotion of exudation	✓				✓			✓			✓
Poss. surgical debridement (then new wound classification), rehydration, promotion of exudation		✓			✓			✓			
Poss. surgical debridement (then new wound classification), exudation management, promotion of autolysis, protection of wound edges		✓				✓		✓			✓
Poss. surgical debridement (then new wound classification), rehydration, promotion of autolysis		✓			✓			✓	✓	✓	
Poss. surgical debridement (then new wound classification), exudation management, promotion of autolysis, protection of wound edges		✓				✓		✓	✓	✓	✓
Rehydration, protection against contamination, ensures nondisturbance of the wound			✓		✓	✓		✓	✓	✓	
Exudate management, protection against contamination, ensures nondisturbance of the wound			✓			✓		✓			✓
Rehydration, protection of wound edges and newly formed epithelium tissue, ensures nondisturbance of the wound			✓					✓	✓	✓	
Exudate management, protection of wound edges and newly formed epithelium tissue, ensures non-disturbance of the wound			✓					✓	✓	✓	✓
Protection of newly formed tissue, protection against mechanical impact								✓			
Infection control, rehydration, odor control, debridement (if needed)				✓	✓		✓				✓
Infection control, exudation management, odor control, debridement (if needed)				✓	✓		✓				✓

Askina® Foam

Highly absorbing foam dressing for heavy exudation



Wound dressings

- Fast absorption, even of thick wound exudate
- Excellent exudate management
- High level of comfort
- Retains shape under compression bandages
- Cost efficient: no need for frequent dressing changes
- Two-layered non-adherent foam dressing

Practical tips on application

Askina® Foam is ideal for use under compression bandages – the retention capacity of the foam is not affected.

In wounds of 1-2 cm depth, wound filler can be dispensed with as the foam structure ideally conforms to the wound bed, providing good contact to the wound base.

Indications

- Surgical wounds
- Traumatic wounds, e.g. cuts and abrasions

Article number	Description	Size	Pcs/Pack
7240703	Askina® Foam	5 x 7 cm	3 pcs
7240710	Askina® Foam	5 x 7 cm	10 pcs
7241003	Askina® Foam	10 x 10 cm	3 pcs
7241010	Askina® Foam	10 x 10 cm	10 pcs
7241210	Askina® Foam	10 x 20 cm	10 pcs
7242005	Askina® Foam	20 x 20 cm	5 pcs

